



揮発性物質用電気炉

VMF

carbolite.com/jp/products/coal-coke-test-equipment/vmf-volatile-matter-furnaces/


VMF

揮発性物質用電気炉

VMFシリーズは、石炭の揮発性試験用にデザインされています。

VMF 10/6


石炭の揮発性試験用電気炉。

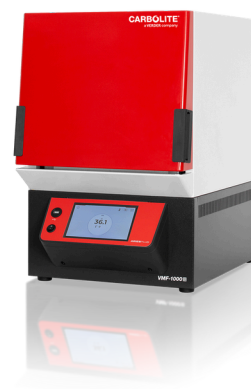

VMF

標準仕様

- VMF / ASTM は、ASTM D3175-11試験法にデザインされた製品です。
- VMF/ASTM、るつぼ、蓋、ホルダー付き
- VMF 10/6 ソフトクロー징、垂直リフトドア
- シースレス熱電対用キャリブレーションポート

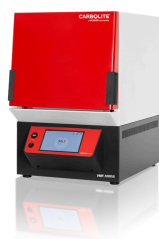

VMF 10/6



VMF 10/6

VMF

オプション (注文時対応)

- アクセサリーとして、坩堝と坩堝用蓋をご用意しています。
- 4個充填坩堝用ラック・取手をご用意しています。(※VMF 10/6のみ)
- 加熱防止機能
- 試験工程において、デジタル通信または nanodac™チャートレコーダー を使用してテスト結果をCSV出力することが可能です。(※VMF 10/6のみ)
- デジタル指示計付き校正済み熱電対プローブ（オプション）（VMF 10/6用）

VMF

るつぼ・ホルダー



オプションの4つのるつぼスタンド
+4つの蓋付きるつぼ (VMF 10/6)



オプションの蓋付きるつぼ (VMF 10/6)



耐腐食性、耐酸化性のインコネル
製るつぼと蓋 (VMF/ASTM) - 本体
に付属



ワイヤーるつぼホルダー

標準的なワイヤーるつぼホルダーは、るつぼと蓋 (VFM/ASTM) と一緒に表示されています - 炉に付属しています。

VMF

製品仕様

VMF 10/6

最高温度 (°C)	1000
連続運転温度 (°C)	950
タイプ	前面投入
昇温時間 (分)	20 to 900°C
寸法: 内部 H × W × D (mm)	100 × 210 × 260
寸法: 外形 H × W × D (mm)	670 × 435 × 608
最大電力 (W)	3000
熱電対タイプ	K
重量 (kg)	47
電源	220 – 240V, 50–60Hz, 単相, 13A

VMF/ASTM

最高温度 (°C)	1000
連続運転温度 (°C)	1000
タイプ	上面投入
昇温時間 (分)	記載なし
寸法: 内部 (mm)	50 × 100
寸法: 外形 H × W × D (mm)	330 × 410 × 300
最大電力 (W)	950
熱電対タイプ	N
重量 (kg)	9
電源	220 – 240V, 50 – 60Hz, 単相 + 接地線, 4A

VMF

資料ダウンロード

製品データシート VMF (日本語) → carbolite.com/files/451168/zhi-pindetashito-vmf.pdf

COAL, COKE & IRON ORE TESTING FURNACES (英語) → carbolite.com/files/17033/coal-coke-iron-ore-testing-furnaces.pdf

VMF

お問い合わせ →

carbolite.com/jp/contact-services/contact-form/?prefilltext=contact_form.request&product=C-2-ASF-B

見積依頼 →

carbolite.com/jp/contact-services/contact-form/?prefilltext=contact_form.quote_request&product=C-2-ASF-B



CARBOLITE®
a **VERDER** company

carbolite.com/jp/products/coal-coke-test-equipment/vmf-volatile-matter-furnaces/

© Carbolite Gero – carbolite.com / Content may be subject to modifications or corrections